

User Manual

Shower Chair with Arms and Back

Model 796A



The user manual contains important safety information. Please read and keep it for future reference.

IMPORTANT SAFEGUARDS



PLEASE READ ALL SAFETY INFORMATION CAREFULLY BEFORE USE.



DO NOT use this product for individuals with severe mental, cognitive, physical, or visual impairments.



DO NOT leave an infected user unattended.



DO NOT stand on the shower chair.



DO NOT use the shower chair for purposes other than intended.



DO NOT place this product on uneven surfaces.



Always ensure that the Legs are properly adjusted to ensure the chair is level and stable in all directions.



Use Armrests for assistance only. **DO NOT** attempt to support the entire body weight with the Armrests.



When using the Backrest, ensure the shower chair remains balanced. Excessive twisting or rotating while sitting on the chair may cause imbalance and product damaged.

*Note

Failure to follow the above noted safety instructions can lead to severe injuries

PARTS LIST




seat panel x1



backrest x1



backrest pad x1

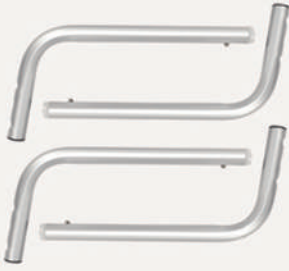

support leg tube x4



armrest x2




seat cushion x1



arc-washer x12



backrest tube x2

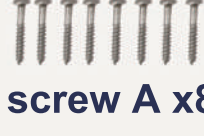

screw A x8



hex-agonal screwdriver x1



storage bagx1



screw B x4



adjustable leg tube x4

1


Step 3

Tighten the leg tubes using screw A, with the arc-washer already attached



Step 4

Seat leg frame installation is complete

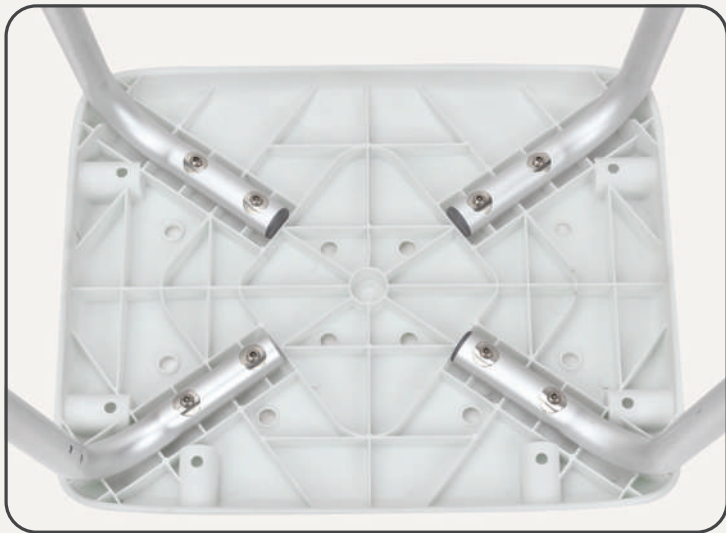

3

Step 7

Insert the armrests into the armrest holes on both sides. Press the pinball with your hands and ensure it pop up through the pinhole



Step 8

Insert cushions as shown in the picture.(Note: The cushion installation is optional and can be skipped if not required.)



5

Step 1

Fix the backrest using screw B with the arc-washer already in place. Tighten it using hexagonal screwdriver



Step 2

Make sure the support legs align with the grooves, then push them in by hand



2


Step 5

Insert the adjustable leg tube into the leg support tube



Step 6

Insert the backrest into the backrest holes. Press the pinball with your hands and ensure it pop up through the pinhole

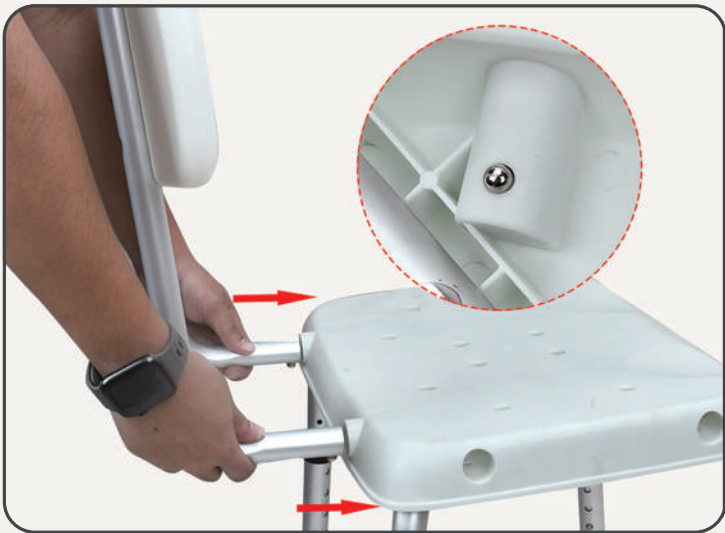

4

Step 9

Install the bag as shown in the picture



MAINTENANCE

1. To better protect the product, clean it with mild soap and warm water after use. Then Wipe it dry and store it in a dry place.

2.Please check whether screws and components are loose regularly, If they are loose, tighten them in time.

6

WARRANTY

Thanks for choosing Lemagalu!

This product is covered with Lemagalu's warranty for 12 months from the date of original purchase.

30-Day Limited Warranty

To provide customers with quality shopping and user experience, we guarantee to replace or make a full refund within the 30-day period from the date of purchase.

30-Day to 1-Year Limited Warranty

After the 30-day warranty has expired, Lemagalu activates an extra 1-year warranty that covers product replacement from the date of purchase.

30-Day to 1-Year Warranty Limitations

The warranty does not apply which damage is caused by other factors, including without limitation:

- (a) Normal wear and tear;
- (b) Misuse, dismantling, accident, or failure to follow operating instructions.
- (c) Exposure to liquid or infiltration of foreign particles;
- (d) Servicing or modifications of the product other than by Lemagalu;
- (e) Commercial or non-household use.

We can only provide after-sales service for products that are sold by Lemagalu or Lemagalu authorized retailers and distributors.

* Please have your order ID ready before contacting Customer Support

7

Need Help? Contact Us!

If you have any questions or concerns about your new product please contact our Customer Support Team.

Lemagalu@outlook.com



How to contact sellers on Amazon ?



Have a question for Lemagalu?

Ask a question

What can I help you with ?

An order I have placed

An item for sale

8